



IMPORTERS / EXPORTERS & SUPPLIERS :-

SUNSHINE INSTRUMENTS

Deals in :- Industrial / R&D / Scientific / Analytical / Diagnostics / Accessories

83, Chinny's Chamber, Dr. Nanjappa Road, Coimbatore - 641 018. (T.N)

Ph : +91 422 2393357, +91 93600 23357 E-mail : instrumentsindia2012@gmail.com

Website : www.lutrontestinginstruments.com



JAYANT SIEVES SHAKERS

CAT No. J - 91 TABLE MODEL

Table Model : Analytical Sieve Shaker are used in research & development, Quality Control of raw material, as well as in Production Monitoring.

The Heavy Cast Iron Unit. The Motor actuates a eccentric weight. This producing uniform motion an sieves table; which is equality balanced mounted on helical springs. The vibration on Sieve Table can be control by speed regulator & imported Mechanical Timer to set required time for sieving analysis.


Specifications :

- | Variable Amplitude by Regulator.
- | On / Off Switch & Pilot Lamp.
- | Mechanical Timer 10 to 60 minutes for Auto - Cut off of set time for Sieve Analysis.
- | L/R required for Sieve Analysis.

Capacity : Maximum 8" x 2" ht Jayant Test Sieves of 6 different Mesh size plus Lid / Receiver.

Weight : Unit : 35 kg Shipping : 45 kg

NOTE : Rates without Test Sieve & Lid / Receiver Wooden Case Packing

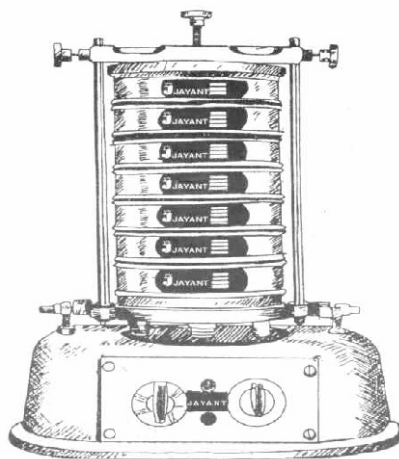

SUNSHINE INSTRUMENTS