



IMPORTERS / EXPORTERS & SUPPLIERS :-

SUNSHINE INSTRUMENTS

Deals in :- Industrial / R&D / Scientific / Analytical / Diagnostics / Accessories

83, Chinny's Chamber, Dr. Nanjappa Road, Coimbatore - 641 018. (T.N)

Ph : +91 422 2393357, +91 93600 23357 E-mail : instrumentsindia2012@gmail.com

Website : www.lutrontestinginstruments.com



CAT NO J - 94 NO-TAP

RO-TAP SIEVE SHAKER : Floor Mode, Heavy Casting Iron Body: No Foundation required.

The laboratory standard in Test Sieves Shaker JAYANT Ro-Tap has a unique two dimensional operation; a horizontal circular motion & vertical tapping.

This unique action allows material particles to stratify & "Seek" Critical Opening in the sieve analysis.

This result is the most accurate & consistent particle analysis testing.

Capacity : To carry maximum 8" DIA Y 2" HT Jayant Test Sieves of 7 different Mesh Size plus Lid / Receiver. Lid / Receiver required with Shaker.

Specifications :

1 Circular (Oscillation) per minute = 300 to 350

Oscillation Displacement : 1" to .75"

Taps per minutes : 180 to 200


Hammer weight : 9.5 kg

Crompton 1/4 HP Motor 230 V, 50 Hz, 180 Watts fitted inside body.

Timer Box : Separated unit built in with : 3 Pin Socket, On / Off Switch, Pilot Lamp & 3 Core Cable with 3 Pin Plug of 15 Amps & Mechanical Timer 10 to 60 Minutes minimum 10 Minutes.

Weight : Unit : 115 kg, Shipping : 155 kg

Note : Rates without Test Sieves & Lid / Receiver Wooden case with Nut / Bolt packing.

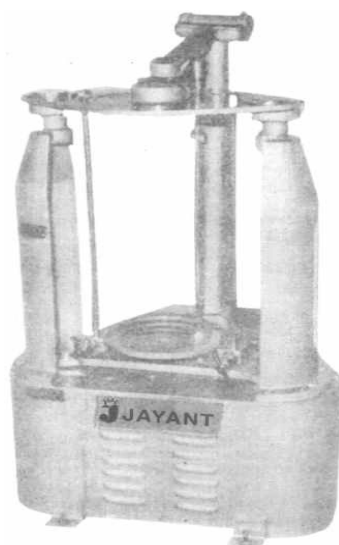

SUNSHINE INSTRUMENTS